




ELECTRIC MASSAGE CHAIR

MASCHR-HOR-B-UP-GYBK

USER MANUAL

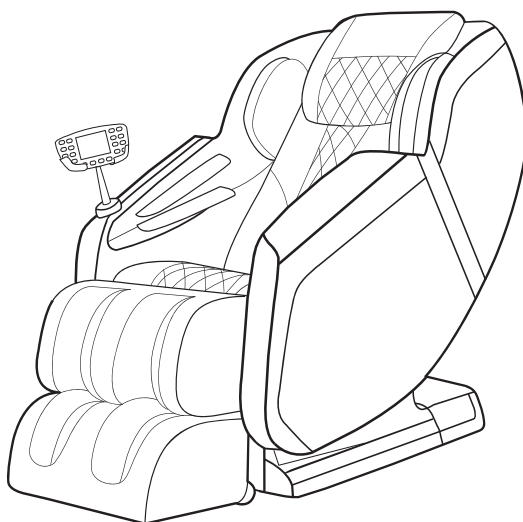


Please read this user manual carefully before using the product.
Retain this manual for future reference.

IMPORTANT SAFETY INSTRUCTIONS

Operating environment

- Do not use or store the massage chair in a damp environment such as a bathroom, sauna or spa to prevent mould, electric shock fire or mechanical failure.
- Do not use the massage chair outdoors.
- Do not use the massage chair under dusty, smelly or anoxic conditions.
- Do not place the massage chair near a heater stove or under direct sunlight where the heat radiated may cause inner components to overheat and cause a fire.
- Ensure that the massage chair is standing on a flat and level ground.
- Ensure that there is an earthed connection for the massage chair to prevent electric shocks.
- Do not use the massage chair while smoking.

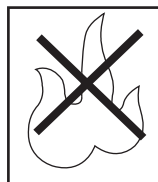

Please consult doctor before using chair if you have the following conditions:

- Suffering from a malignant tumour
- Serious heart disease
- Serious skin disease
- Serious boned hyperplasia
- Resting under doctor's supervision
- Bad back from accident or illness
- Need of peaceful surroundings
- Arthritis
- Implanted pacemaker
- Bowel ailments (eg, gastritis, enteritis, diarrhoea, hepatitis)
- Serious blood circulation problem
- Consciousness problem
- Hyper pyretic

WARNING

Children under 18 or physically or mentally-challenged persons are not allowed to operate the massage chair to prevent electric shock, injuries or mechanical malfunctions.

For safety reasons and to highlight the potential dangers of improper use of the product, the instructions below should be closely adhered to.

Non-professional technicians are not allowed to dismantle or repair the chair to avoid fire, electric shock or serious injuries.	
Do not wet the switch when cleaning so as to avoid an electric shock or short circuit.	
During storage, cleaning or maintenance, ensure that the power plug is disconnected from the socket. Do not detach the plug with wet hands. Make sure that power cord is detached when moving the chair.	
Do not place the chair on an unstable surface to prevent electric shock or chair overturning.	
Clean power plug periodically to prevent fire due to dust.	
Do not use the chair at unsafe places or for other purposes. This may cause product damage or accidents.	
Keep the power cord free of entanglement at all times and do not place heavy items on it.	
The chair has a heated surface. Ensure that persons who are insensitive to heat to be careful when using the chair.	

Do not modify the structure or change components by unauthorised persons. It may cause an electric shock, short circuit or fire.	
Do not store the chair near fire or flammable substances. It may cause product damage and defects.	
Should any problem or abnormality (such as burning smell) occur, stop using the chair immediately and detach the power plug from the socket. Operating under such circumstances may cause a fire or electric shock. Contact your supplier or authorised agent for repairs and advice.	
Only use the rated voltage specified for the chair.	
Stop using the chair immediately if the power cord or power plug is damaged or loose.	
Keep children away from the chair when it is in operation to avoid any accident or injury.	
Do not place the chair under direct sunlight at all times.	
If the power cord is damaged, it must be replaced by the manufacturer, supplier or authorised agent.	

WARNING

Supervision and guidance are necessary when children, physically or mentally-challenged persons, and persons inexperienced or lack of knowledge with the use of massage chairs intend to use the massage chair. Unless there is someone who is able to be responsible for their safety, the massage chair should be unavailable to them. Please ensure that children do not play with the chair. The chair has a heated function and must be used with caution with anyone who is insensitive to heat.

SPECIFICATIONS


Product name	Electric Massage Chair
Model	MASCHR-HOR-B-UP-GYBK
Power supply	240V 50Hz
Rated power	90W
Rated working time	15min




WARNING

Use only the rated supply of 240V to avoid damage to electrical and electronic components of the chair.

MASSAGE CHAIR PARTS

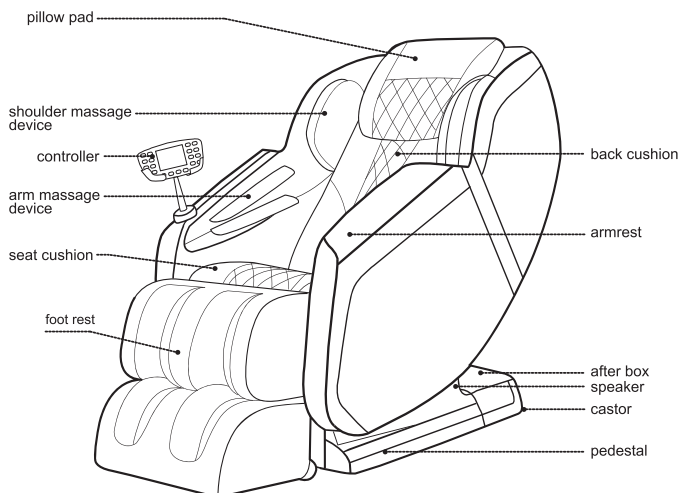

REMOTE CONTROL FUNCTIONS

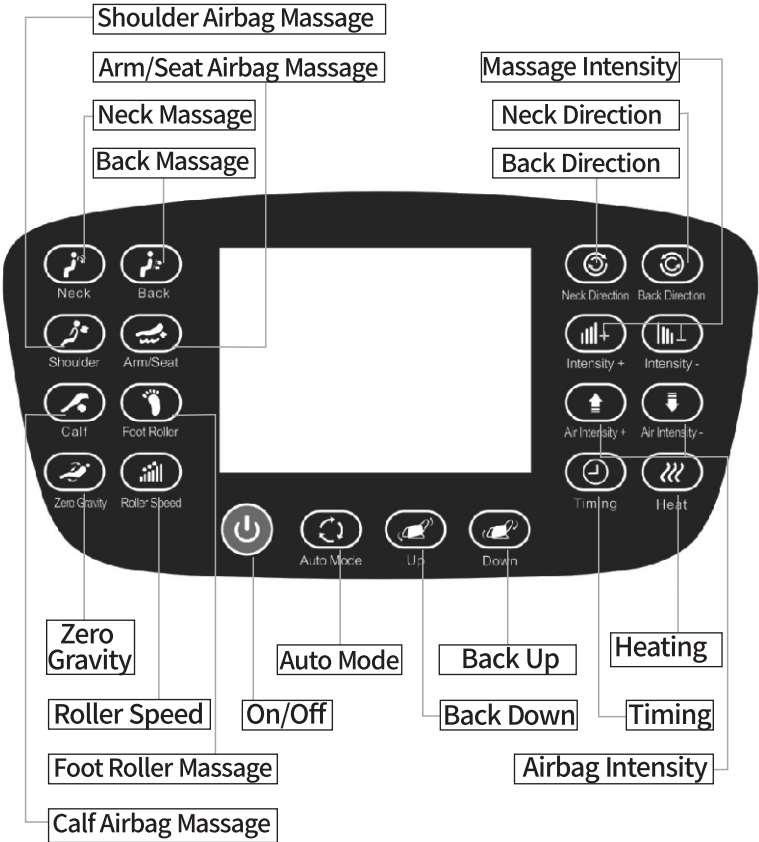


Warning: Only rated supply 240V~ must be used to avoid damage to electronic elements.

OPERATIONS GUIDE

1. Turn on the power and turn the switch to the “I” position.

- ▶ Insert the power plug into power socket, and push the power switch to “I” position. The chair is in resetting condition, remote power key indicator will flash quickly. After reset, the chair goes to standby mode.

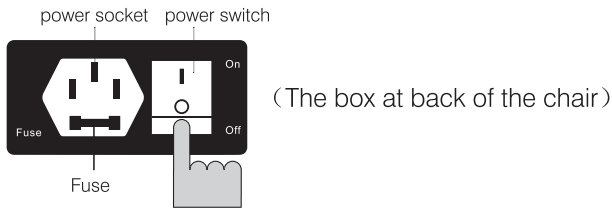

Warning: Do not insert or pull out the power plug with wet hands, otherwise an electric shock would happen.

2. Sit on massage chair and press the “ON” button to start.

Remote Control Usage Guide



•On/Off

Press this button to start/stop the massage system.



•Auto Mode

Press this button to switch an auto program each time, there are 20 auto programs for choice



•Neck Massage

Press this button to start/stop the neck massage



•Back Massage

Press this button to start/stop the back massage



•Neck Direction

Press this button to choose the neck massage direction



•Back Direction

Press this button to choose the back massage direction

- 

•Roller Speed
Press this button to switch foot roller speed each time ,3 levels of massage speed for choice
- 

•Heating Function
Press this button to start/stop the heating in the back and calf
- 

•Foot Roller Massage
Press this button to start/stop the foot roller massage
- 

•Shoulder Airbag Massage
Press this button to start/stop shoulder airbag massage
- 

•Arm/Seat Airbag Massage
Press this button to start/stop arm or seat airbag massage and seat roller massage.
- 

•Calf Airbag Massage
Press this button to start/stop calf airbag massage
- 

•Back Up
Press this button, the backrest will be moved up to the original position,release this button, it will stop.
- 

• Back Down
Press this button, the backrest will be moved down to the lowest position,release this button, it will stop.
- 

• Massage Intensity
Under manual mode,Press“  ”button to increase the massage speed,Press“  ”button to decrease the massage speed.
6 levels of intensity for choice.
- 
- 

• Airbag Intensity
Open airbag function,Press“  ”button to increase the airbag intensity,Press“  ”button to decrease the airbag intensity.
6 levels of intensity for choice.
- 
- 

• Timing
Press this button to set the timing of 15-20-25-30 minute
- 

• Zero Gravity
Press this button, the chair will go to zero-gravity position.
Press it again, the chair will go back to original position.

Mode Switching

Auto Mode: In auto mode, there are 20 auto programs for choice.

Manual Mode: In manual mode, choose any of manual function on the remote control to have a specific massage.

Bluetooth Connection

Bluetooth function is turned on when the machine turned on.

Turn on your mobile Bluetooth function and search the massage Bluetooth: "livemor" and connect to play the music.

3 Restore massage chair.

Press "  " button after using.

Warning

Do not stand in the leg cushion before it restores to original position to prevent any damage to leg cushion.

 **Warning:** Supervision and guidance is mandatory for those disabled mentally and physically. Please advise with your doctor before using the chair if you have ANY medical conditions. Children under 18 should not use the chair and consult with a doctor for approval if attempted. Do not operate the chair if there are any children or pets near by. Serious injuries may occur if these instructions are not followed.

SAFETY INSTRUCTIONS

- Should any abnormality or problem occurs with the chair, switch off the chair and detach the power cable from the power socket.
- When the power cord is damage, it must be replaced by one from the manufacturer or agent.
- Do not attempt to dismantle, modify or change components by unauthorised persons to avoid any electric shock, fire or injury.

CARE AND MAINTENANCE

Storage

- Do not dismantle the back of the chair. Electric shocks may occur when touching inner components of the chair.
- Do not dismantle the chair by unauthorised personnel to prevent electric shocks, fire and injury.
- Ensure that the chair is switch off and the power plug is detached from the socket when cleaning the chair to avoid any electric shock or hazard.
- Keep the massage chair clean and free of dust. Do not place the chair in a damp environment.
- If the chair is not in use for long periods, please cover it with a cloth to protect it from dust.
- Do not place the chair under direct sunlight.

Maintenance

- Use a neutral detergent only to clean the back cushion, head cushion and leg cushion. Let it air dry.
- Use a cloth and detergent to clean the plastic parts, pipe and leather material of the massage chair.
- Do not use organic solvents such as gasoline, amyl acetate or benzene to clean to prevent leather colour fading, ageing or damage.
- Do not spray pesticides on the chair to avoid electric shocks or fire.
- Only clean with a dry cloth to prevent mechanical failure.

Remote controller care

- Do not place the remote controller or massage chair near to a heater or fire source.
- Remove dust on remote controller with a soft cloth or small banister brush.
- Do not use volatile solvents like benzene, gasoline or amyl acetate to clean the remote controller.

TROUBLESHOOTING

Check the massage chair periodically when it is in frequent use for any defects or damage.

If the massage chair has not been used for some time, please ensure that all components are in working order before using it.

PROBLEM	SOLUTION
Failure to start and indicator light is off	1. Check that plug is well connected. 2. Check that switch on chair is set to “I” position.
Sounds during operation	These are normal sounds during operations 1. Sounds from kneading or knocking. 2. Sounds from speed adjustment. 3. Sounds from box in massage mode. 4. Sounds from vibration in back massage. 5. Sounds from air input and output. 6. Sounds from airbag inflating and venting. 7. Sounds from leg adjustments. 8. Sounds from motor when massage chair is loaded.

