

Livemor™

Massage Chair

MASCHR-G-BERYL-GR

An instruction manual



Pictures are for reference only

Living Art

Thank you for choosing our massager products.

This massage chair adopts ergonomic design to massage the head, neck, back, shoulder, arm, waist, hip, leg, foot and other parts, with various massage methods, warm compress, air bag extrusion massage and other functions. It is a luxury high-end massage chair.

Please read this manual carefully before use and pay special attention to safety precautions.

To ensure proper operation and use of this product.

Keep the manual safe for future inquiries.

Catalog

I. Product Characteristics

II. Usage

III. Maintenance and Maintenance

IV. Matters needing attention

V. Safety Protection

VI. Warm Tips

VII. Technical parameter

VIII. Circuit schematic diagram

IX. Storage and Maintenance

X. After-sales Service Guarantee Card

XI. Analysis and Removal of Common Faults

XII. Analysis and Removal of Common Faults

Failure phenomenon	Possible cause	Processing Method
Sound (working sound) *inflatable sound of airbags *bleed sound of airbags *changing-over sound on solenoid valve *air pump working sound	The sound produced by the product at work Sound of mechanical operation	This sound is normal No need to deal with it
The whole machine is not working at all	*The power plug is not inserted into the socket *Socket switch not open *The massage button is not pressed	Check the power plug and socket, ensure that the plug is in good contact with the socket, press the massage button.
The massage stopped suddenly in the process	This product has just completed 15 minutes of rated working time and automatically switched to standby mode.	If you would like to massage again, please do so again. Press the  key to turn on massage.
	Overheating protection for motors exceeding rated service time	Turn off , Cut the power, use the unit again after the motor is cool down

*The company reserves the right to improve products and instructions.

If alteration is made without notice, the appearance of the product shall be in accordance with material object.

*This volume has been carefully checked, but it does not exclude the possibility of misprints or misunderstandings in the content.

We are sorry to inform you that we reserve the right of interpretation.

 Improper disposal of waste causes environmental pollution.

X. After-sales Service Guarantee Card

Warranty Card

With this card warranty for one year, users will keep it properly.

User Union

Buyer ID: _____
 Order No.: _____
 Product Name: _____
 Model No.: _____
 Store: _____
 Date: _____ Price: _____
 User name: _____
 Telephone: _____
 Address: _____

Warranty Card

With this card warranty for one year, users will keep it properly.

Company Stubs

Buyer ID: _____
 Order No.: _____
 Product Name: _____
 Model No.: _____
 Store: _____
 Date: _____ Price: _____
 User name: _____
 Telephone: _____
 Address: _____

I. Product Characteristics

- The whole machine has beautiful lines, bright colors, humanized design and green environmental protection materials.
- This product is designed with different groups of movements and many massage heads, which can knead and massage the neck, back, waist and buttocks of humanbody.
- The whole machine is equipped with about many air bags divided into legs, arms, cushions and shoulders to freely display the function of air bag finger pressing and massage.
- It can manually adjust the massage movement of neck, back, waist and buttock. Various massage techniques can be changed freely, and three-level adjustment can be realized.
- The recline angle of the back can be adjusted according to the needs, When the angle of the back is adjusted to the minimum, the zero gravity function can be realized.
- One key heating.
- The massage time can be adjusted according to the needs, with 5 minutes plus and minus steps, and a massage journey can be set for 15 minutes, 20 minutes, 25 minutes and 30 minutes.
- Specially designed buttocks and sole rollers.
- Special hifi Bluetooth audio, massage at the same time, enjoy the fun of music.

Maintenance record:

Date	Fault Condition	Address	Result

II. Usage

1. Installation steps of manual support

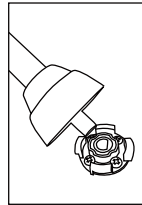


Figure 1

Align and insert the plane of the plug at the bottom of the manual control bracket with the groove plane of the bracket fixing device above the right armrest (as shown in Figure 1), and then insert it to the end. Finally, turn the round knob at the bottom of the bracket clockwise until it is tight.

2. Connecting wire installation method



Figure 2

Connect the pin plug of the power line with the plug of the electrical box (as shown in Figure 2), and connect the other end with the power socket.

IX. Storage and Maintenance

Cleaning

Never use steel wool, abrasive cleansers or corrosive liquid to clean your product.

❶ Be sure to disconnect the power cord from the electrical outlet before cleaning.

❷ Clean the massage chair, remote controller and power cord with a soft dry cloth.



❸ Do not use wet cloth or detergent to clean the product.

Storage

❶ Disconnect the power cord from the electrical outlet.

❷ Please keep the product away from direct sunlight.

❸ Do not stack heavy items on the product.

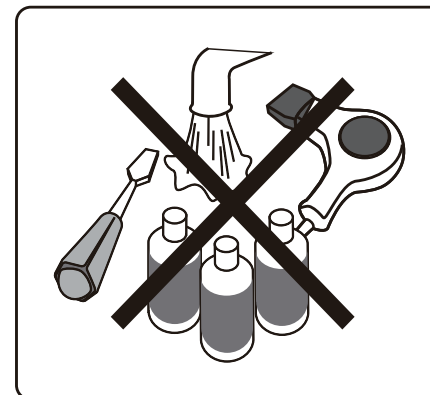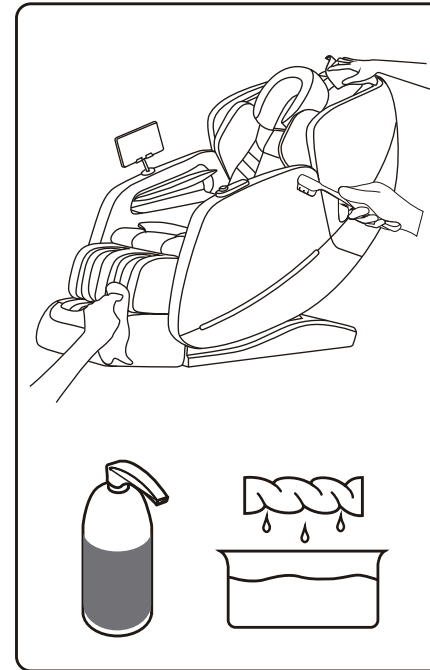


Never operate this product if it has damaged/exposed parts, please contact the authorized dealer for maintenance.

Please disconnect the power cord from the electrical outlet when the massage chair sits idle for extended periods of time.


VII. Technical parameter

Product Name:	Massage chair
Rated voltage:	240V~
Rated frequency:	50Hz
Rated power:	80W
Rated working time:	15mins.
Security level:	Class I

VIII. Circuit schematic diagram

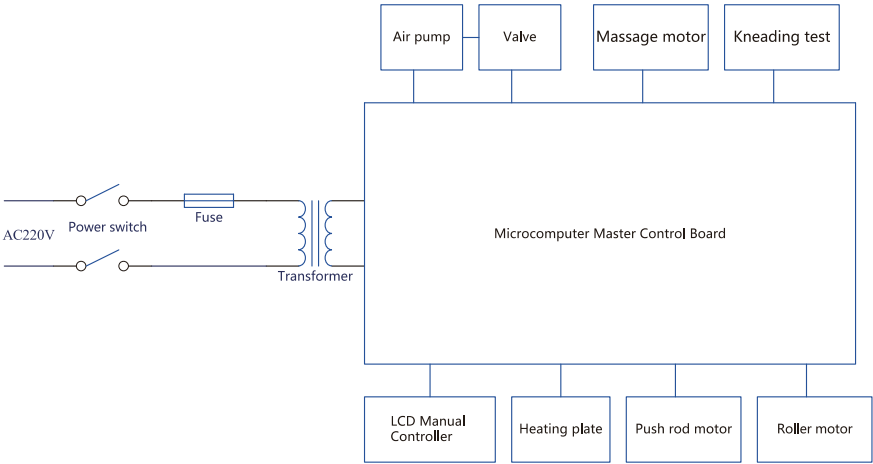

3.Manual operation instructions

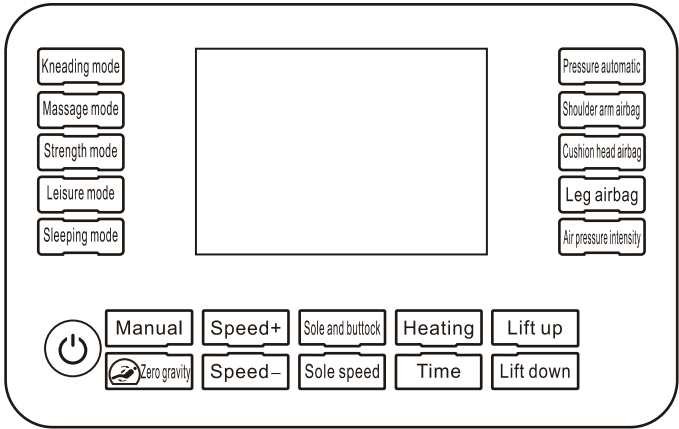

- : Start or close.
- Kneading mode: neck back kneading massage, and then press this key message head is reversed adjustment.
- Massage mode: waist back push to take massage, and then press this key message head is reversing adjustment.
- Strength mode: Fully automatic power mode.
- Leisure mode: Fully automatic leisure mode.
- Sleeping mode: Fully automatic sleep mode.
- Pressure automatic: Turn on or off the whole body airbag.
- Shoulder arm airbags: Turn shoulder and arm airbags on or off.
- Cushion head airbags: Turn seat cushions, head, leggings, airbags on or off.
- Leg airbags: Turn the leg airbags on or off.
- Air pressure intensity: Adjust the air pressure intensity.
- Manual: Stop all massages.
- Zero gravity: comfortable degrees elevated between the body and legs.
- Speed plus: Massage head speed plus.
- Speed-: Massage head speed-.
- Sole and buttock: Press the sole and buttock rollers for synchronous massage, the second lower hip roller massage, and the third foot roller massage.
- Sole speed: The sole roller speed adjustment.
- Heating: Waist,calf.
- Time: Time setting, 5 minutes for each click, up to 30 minutes.
- Lift up: Long press back to pick up (massage movement, air pump will stop working when lifting).
- Lift down: Press and hold back (massage movement, air pump will stop working when down).

III. Maintenance and Maintenance

1. Please use neutral detergent and soft cloth for daily cleaning, and do not use corrosive detergent.
2. If the power cord is damaged, it must be replaced by the guidance of professional personnel to avoid danger.
3. When not using this product or cleaning, please unplug the power supply to avoid danger.
4. Do not use any sharp objects to cut the product surface.
5. Do not put the product in a humid environment, but in a dry and ventilated place.
6. Do not press or press the product hard.

IV. Matters needing attention

In order to ensure the safety and correct use of the massage chair, please refer to the operation manual in detail before using this product.

To reduce the risk of fire, burns, electric shock, or injury to others, follow these instructions:

- 1) Power supply other than AC 220V-240V shall not be used.
- 2) When no one is using the machine, please be sure to remove the power plug of the massage chair from the socket. When not using the machine or before removing the device / parts, please also remove the power plug from the socket in time.
- 3) In order to avoid fire, electric shock or injury caused by overheating, please do not put the product under pillows, blankets and other covers.
- 4) If the product is used for children or the disabled, adult guardians are required to supervise and guide it.
- 5) Please be sure to use the product according to the usage described in the user manual. Do not use parts that are not recommended by the manufacturer.
- 6) Do not use this product in case of any of the following conditions: the power cord or plug is damaged, the operation is not normal, if it falls to the ground or other damages, or falls into water, please send the product back to the service center for inspection and maintenance.
- 7) Do not drag the power cord to move the product or use the power cord as a handle.
- 8) Keep a safe distance between the power cord and the heating surface.
- 9) Do not operate the product when the air outlet is blocked. Please keep the air outlet unobstructed to avoid being blocked by cotton chips, hair and other things.
- 10) Do not insert or drop objects into any openings of the product.
- 11) Do not use this product outdoors.
- 12) please do not operate the product in the place of spray product or oxygen filled place.
- 13) Cut off the power supply, please turn off all operation controls, and then unplug the power plug from the socket.
- 14) Do not cover the massage chair so as not to cause malfunction or fire due to excessive temperature.

V. Safety Protection

- Do not use the massage chair in excessively humid or dusty environment to avoid electric shock or mechanical failure.
- Do not use the massage chair when the room temperature exceeds 40 °C.
- Do not place the massage chair near the electric heater, stove or direct sunlight.
- The massage chair is designed for personal use and is not suitable for any commercial use.

Careful

- Massage chair is only for indoor use.
- Use massage chair at low room temperature, do not suddenly adjust the temperature. It is recommended to gradually adjust the temperature to the normal range before use.
- If the massage chair is moved from a low temperature storage place to a warm environment for use, it may affect the normal operation of the machine due to the accumulation of moisture in the parts, or even cause parts failure. It is recommended that the product be placed at normal room temperature for 1 hour before use.

VI. Warm Tips

- The power switch must be turned on before other functions can be selected.
 - Children should use this massage chair under the supervision of adults to avoid equipment damage and children's injury due to improper operation.
 - Please do not let children climb this massage chair to avoid children falling.
 - Use the massage chair to adjust the strength of the chair massage, because the massage force is too small to play its due role in stimulating, too large prone to fatigue, and damage the skin, massage force and times from less to more, from light to heavy.
 - Use the massage chair to relax the body and mind. In addition to thinking at home and abroad, you should be calm and relaxed.
 - Massage chairs are not suitable for patients with heart disease, hypertension, osteoporosis, local skin damage, ulcer, bleeding, tuberculosis, tumor, people who are too hungry, full, drunk or tired.
 - In order to avoid the reduction of massage effect or injury, please do not use the massage chair together with other effective equipment such as electric blanket.
 - Do not use the massage chair within 1 hour after eating to avoid adverse effects.
 - Do not use the massage chair when your hands or other parts of your body are wet. Do not plug or unplug the power plug.
 - Do not insert any articles or inflammables into the massage chair. This may cause electric shock, malfunction or fire.
 - If you feel any discomfort during massage, please stop immediately.
 - If the product is not used for self diagnosis or treatment and needs diagnosis or treatment, it shall be handled by professional medical personnel.
 - Do not use this product if you are in medical condition and are pregnant or in physiological period.
 - Use the massage chair to master the massage time, 15 minutes for each massage, and no more than half an hour for each person (time) (when it reaches half an hour, let the massage chair rest for ten minutes before use). If it is used continuously for a long time, the overheat protection device may be started and the power will be cut off automatically.
 - The surface of the appliance is hot. People who are not sensitive to the heat must pay attention to it when using it.
 - Do not directly massage elbow joint, knee, chest, abdomen, etc.
- Not suitable for massage chair:**
- In case of any of the following conditions, please turn off the switch immediately, cut off the power supply and contact our company.**
- accidentally spilling liquid on the massage chair may cause electric shock.
 - the cloth cover of massage chair is damaged and the internal parts are exposed.
 - in case of abnormal or severe pain during massage, please stop massage immediately and seek medical consultation.
 - please stop using immediately if any mechanical failure or abnormal condition is detected during use.
 - in case of current interruption in use, attention shall be paid to avoid user injury due to sudden power recovery.
 - in case of lightning.