EQUIPMED®



Colours and designs may vary

Bath Lift Chair

User Manual

[Revision 1.0]

READ THIS MANUAL CAREFULLY BEFORE USE – FAILURE TO DO SO MAY RESULT IN INJURY, PROPERTY DAMAGE AND MAY VOID WARRANTY. • KEEP THIS MANUAL FOR FUTURE REFERENCE. • Products covered by this manual may vary in appearance, assembly, inclusions, specifications, description, and packaging.

THE BATTERY MAY HAVE SOME CHARGE ON ARRIVAL, IT IS IMPORTANT HOWEVER THAT THE BATTERY IS CHARGED PRIOR TO USE.

Safety

Safety messages are designed to alert you to possible dangers or hazards that could cause death, injury or equipment or property damage if not understood or followed. Safety messages have the following symbols:



You WILL be KILLED or SERIOUSLY INJURED if you do not follow instructions.



You CAN be KILLED or SERIOUSLY INJURED if you do not follow instructions.



You CAN be INJURED if you do not follow instructions or equipment damage may occur.

It is vital that you read and understand this user manual before using the product, including safety warnings, and any assembly and operating instructions. Keep the manual for future reference.

Safety precautions and recommendations detailed here must be fully understood and followed to reduce the risk of injury, fire, explosion, electrical hazard, and/or property damage.

Safety information presented here is generic in nature – some advice may not be applicable to every product. The term "equipment" refers to the product, be it electrical mains powered, battery powered or combustion engine powered.

- **Before Use** - If you are not familiar with the safe operation/handling of the equipment or are in any way unsure of any aspect of suitability or correct use for your application, you should complete training conducted by a person or organization qualified in safe use and operation of this equipment, including fuel/electrical handling and safety.
- Do NOT operate the equipment in flammable or explosive environments, such as in the presence of flammable liquids, gases, or dust. The equipment may create sparks or heat that may ignite flammable substances.
- Keep clear of moving parts.
- Equipment may be a potential source of electric shock or injury if misused.
- Do NOT operate the equipment if it is damaged, malfunctioning or is in an excessively worn state.
- Do NOT allow others to use the equipment unless they have read this manual and are adequately trained.
- Keep packaging away from children - risk of suffocation! Operators must use the equipment correctly. When using the equipment, consider conditions and pay due care to persons and property.

General Work Area Safety

- Work areas should be clean and well it.
- Do not operate the equipment if bystanders, animals etc are within operating range of the equipment or the general work area.
- If devices are provided for connecting dust extraction / collection facilities, ensure these are connected and used properly. Dust collection can reduce dust-related hazards.

General Personal Safety

- Wear appropriate protective equipment when operating, servicing, or when in the operating area of the equipment to help protect from eye and ear injury, poisoning, burns, cutting and crush injuries. Protective equipment such as safety goggles, respirators, non-slip safety footwear, hard hat, hearing protection etc should be used for appropriate equipment / conditions. Other people nearby should also wear appropriate personal protective equipment. Do not wear loose clothing or jewellery, which can be caught in moving parts. Keep hair and clothing away from the equipment.
- Stay alert and use common sense when operating the equipment. Do not over-reach. Always maintain secure footing and balance.
- Do not use the equipment if tired or under the influence of drugs, alcohol, or medication.
- This equipment is not intended for use by persons with reduced physical, sensory, or mental capabilities.

General Fuel Safety

- Petrol/fuel/gasoline is extremely flammable – keep clear of naked flames or other ignition sources.
- Do not spill fuel. If you spill fuel, wipe it off the equipment immediately – if fuel gets on your clothing, change clothing.
- Do NOT smoke near fuel or when refuelling.
- Always shut off the engine before refuelling.
- Do NOT refuel a hot engine.
- Open the fuel cap carefully to allow any pressure build-up in the tank to release slowly.
- Always refuel in well-ventilated areas.
- Always check for fuel leakage. If fuel leakage is found, do not start, or run the engine until all leaks are fixed.

General Carbon-Monoxide Safety

- Using a combustion engine indoors **CAN KILL IN MINUTES**. Engine exhaust contains carbon-monoxide – a poison you cannot smell or see.
- Use combustion engines OUTSIDE only, and far away from windows, doors, and vents.

General Equipment Use and Care

- Keep packaging away from children and pets, potential risk of suffocation or choking.
- The equipment is designed for domestic use only.
- Handle the equipment safely and carefully.
- Before use, inspect the equipment for misalignment or binding of moving parts, loose components, damage, or any other condition that may affect its operation. If damaged, have the equipment repaired by an authorised service centre or technician before use.
- Prevent unintentional starting of the equipment - ensure equipment and power switches are in the OFF position before connecting or moving equipment. Do not carry equipment with hands or fingers touching any controls. Remove any tools or other items that are not a part of the equipment from it before starting or switching on.
- Do not force the equipment. Use the correct equipment for your application. Equipment will perform better and be safer when used within its design and usage parameters.
- Use the equipment and accessories etc. in accordance with these instructions, considering working conditions and the work to be performed. Using the equipment for operations different from those intended could result in hazardous situations.
- Always keep equipment components (engines, hoses, handles, controls, frames, housings, guards etc) and accessories (cutting tools, nozzles, bits etc) properly maintained. Keep the equipment clean and, where applicable, properly lubricated.
- Store the equipment out of reach of children or untrained persons. To avoid burns or fire hazards, let the equipment cool completely before transporting or storing. Never place or store the equipment near flammable materials, combustible gases, or liquids etc.
- The equipment is not weather-proof, and should not be stored in direct sunlight, at high ambient temperatures or locations that are damp or humid.
- Do not clean equipment with solvents, flammable liquids, or harsh abrasives.
- For specific equipment safety use and care, see Equipment Safety.

© Equipmed

General Electrical Safety

- Inspect electrical equipment, extension cords, power bars, and electrical fittings for damage or wear before each use. Repair or replace damaged equipment immediately.
- Ensure all power sources conform to equipment voltage requirements and are disconnected before connecting or disconnecting equipment.
- When wiring electrically powered equipment, follow all electrical and safety codes.
- Always use a Residual Current Device (RCD) (GFCI) / Earth Leakage Circuit Breaker / Safety Switch when operating this item (electrical items only).
- High voltage / high current power lines may be present. Use extreme caution to avoid contact or interference with power lines. Electrical shock can be fatal.

General Electrical Safety

- Electrically grounded equipment must have an approved cord and plug and be connected to a grounded electrical outlet.
- Do NOT bypass the ON/OFF switch and operate equipment by connecting and disconnecting the electrical cord.
- Do NOT use equipment that has exposed wiring, damaged switches, covers or guards.
- Do NOT use electrical equipment in wet conditions or in damp locations.
- Do NOT use electrical cords to lift, move or carry equipment.
- Do NOT coil or knot electrical cords, and ensure electrical cords are not trip hazards.

General Service Information

- The equipment must be serviced or repaired at authorised service centres by qualified personnel only.
- Replacement parts must be original equipment manufacturer (OEM) to ensure equipment safety is maintained.
- Do NOT attempt any maintenance or repair work not described in this manual.
- After use, the equipment and components may still be hot – allow the equipment to cool and disconnect spark plugs and/or electrical power sources and/or batteries from it before adjusting, changing accessories, or performing repair or maintenance.
- Do NOT adjust while the equipment is running.
- Perform service-related activities in suitable conditions, such as a workshop.
- Replace worn, damaged or missing warning/safety labels immediately.

Child Choking Hazard!



- **WARNING: CHOKING HAZARD - CHILDREN UNDER 8 YRS. CAN CHOKE OR SUFFOCATE ON SMALL PARTS. ADULT SUPERVISION REQUIRED.**
- Carefully inspect anything a baby or child could gain access to, including this product & packaging.
- Dangers include anything of a size that could become stuck or lodged in a baby or child's airway.
- Dangers include things like raw, hard fruit and vegetable pieces, large pieces of meat, bones or sausage skins, popcorn, nuts, hard lollies and corn chips, small magnets and batteries, coins, beads, marbles and small uninflated balloons, broken toys and smaller toys, buttons, keys, and coins. But there are also many more.
- Pay particular attention to accessories, nuts, bolts, screws, washers, caps, covers, loose parts, parts that could be removed or break off, pieces of packaging, staples, tape etc.
- In the case of any doubt, secure the item of concern in an appropriate manner, or completely remove the possibility of access by a baby or child.


Before Operation

- Before first use, read these instructions carefully to become familiar with the bath lift and its operation.
- Keep these instructions for anyone who may use or handle the bath lift.
- Before use, assess whether the bath lift is suitable for the user. Consider:
  - user weight
  - body size
  - cognitive ability
  - joint mobility
- Do not allow unsuitable users to use the bath lift.
- If the user has a disability:
  - an assistant must be present
  - an alarm device should be within reach
- Use the bath lift only in a standard bath that meets the required conditions.
- Ask an assistant or qualified professional to confirm that the bathroom environment is suitable.
- Before each use:
  - check that the suction cups and seat surface are clean
  - remove any dirt, stains or residue
  - check that all parts are fitted, complete and undamaged
  - perform a trial operation
  - confirm that the emergency stop button works correctly
- If using bath oils, bath salts or other toiletries, ensure they do not build up or affect the frame or sliding function of the bath lift.
- The user should sit down slowly and calmly to avoid injury or damage to the bath lift.
- Operate the bath lift only with the supplied hand control.
- Install the bath lift according to the installation instructions.

© Equipmed

Prohibited Actions

- Do not exceed the maximum load.
- Do not stand or kneel on the bath lift.
- During use, do not sit on the opposite side of the seat.
- Do not use the bath lift on steeply sloped surfaces over 5° or on stepped surfaces.
- Do not use the bath lift if any parts are missing.
- Do not use the bath lift for any purpose other than bathing.
- Do not leave children unattended near or with the bath lift.
- During use or maintenance, keep hands, feet, other body parts and objects away from moving parts, including:
  - the hinge between the back plate and base plate
  - the hand control hook on the back plate
  - the lift bracket on the base plate

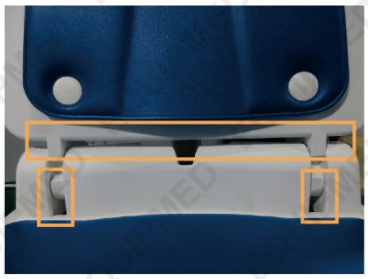

①Hinge between back plate and base plate



②Back plate——hook of controller



③Base plate——at the base plate lift bracket

- Do not place small items, such as toothbrushes or razors, near the cross structure.
- Do not place the hand control near heat sources.
- Do not modify this device.

Important Safety and Battery Information

- Use only the supplied charger to charge the hand control.
- The battery may reach its optimal capacity after several charging cycles.
- During charging, keep the charger and hand control away from:
  - water
  - heat sources
  - damp or wet areas
- The hand control has three-colour indicator lights:
  - Green light: sufficient battery power.
  - Blue light: low battery power. Charge the hand control.
  - Red light: critically low battery power. Charge the hand control immediately to prevent sudden power failure during use.
- During normal charging:
  - the charger indicator light is red while charging
  - the charger indicator light turns green when fully charged
  - charging takes approximately 2.5 hours
- The charger is designed only for charging the hand control. Do not use it for any other purpose.
- The battery is built into the hand control. Do not disassemble the hand control.
- Before transporting the bath lift:
  - disconnect the cables
  - ensure the device is not charging
  - secure all parts to prevent damage or injury caused by falling or impact

© Equipmed

- Before cleaning the hand control or charger:
  - disconnect all plugs
  - clean the charger only with a dry cloth
- Although the hand control and charger are designed to be reliable, they may be damaged if dropped.
- If there are signs of damage on the hand control or charger, stop using them and contact your dealer or the manufacturer.
- Sit in the centre of the bath lift seat.
- Lean back against the backrest.
- Stretch your legs forward.
- Place one hand on the side panel or your lap, and operate the hand control with your other hand.
- The bath lift has passed EMC testing. To reduce the risk of electromagnetic interference, avoid operating the bath lift at the same time as other electrical appliances or a mobile phone in the bathroom.
- If the bath lift operates abnormally, stop using it and contact your dealer for assistance.

Important Information About the Hand Control Battery

- The hand control contains a rechargeable lithium-ion battery.
- Charge the hand control before first use.
- Lithium-ion batteries must not exceed 30% of their nominal capacity during transport. Charge the hand control for at least 12 hours before first use.
- When storing the hand control, keep the battery level at 40–60%.
- For long-term storage:
  - disconnect the hand control from the bath lift
  - store it in a cool, dry and clean environment
  - keep it away from dust and corrosive gases
- The recommended storage conditions are shown below.

Storage time	Storage environment
1 year	10–25 °C, 40–70% relative humidity, non-condensing
3 months	0–35 °C, 30–75% relative humidity, non-condensing
1 month	-10–50 °C, 20–80% relative humidity, non-condensing

- Charge and discharge the battery within 6 months of receiving the product.
- During storage, charge and discharge the battery at least once every 6 months.
- Rechargeable batteries are consumable items with a limited service life. Battery performance will gradually decrease over time.
- Normal battery wear and tear is not covered under returns or claims.
- The manufacturer provides a one-year warranty for the battery.

Safety Symbols

The product may have safety warning labels attached to it, explained below. Understand the symbols on your product and their meanings. If any stickers become unreadable, unattached etc., replace them.



Flammable Material Hazard
Flammable liquids, gases, or substances etc may present. Avoid ignition sources and open flames. Danger of fire.



Read User Manual
Read and fully understand product safety warnings, operation, procedures etc before using the product.



Use Hand Protection
Wear appropriate hand protection and take due care as the product or use of the product may present hand hazards.



Carbon-Monoxide Hazard
Do not use the product in confined areas or without adequate ventilation. Carbon-monoxide poisoning can be fatal.



Electrocution / Electrical Shock Hazard
High voltage or high current electricity may be present or required by the product. Take due care when handling electrical products, cables, plugs and leads. Electrical shock can be fatal.



Toxic Fumes / Dust Hazard
Using the product or by-products from use may produce fumes, smoke or particles that could be harmful if inhaled. Wear appropriate breathing protection and have adequate ventilation.



Explosive Material Hazard
Combustible liquids, gases, or substances etc may be present. Avoid ignition sources and open flames. Danger of explosion.



Cutting / Amputation Hazard
The product may have blades, edges or mechanical devices that can cause severe cut injury to fingers, limbs etc. Take due care when handling and using the product.



Crush Hazard
The product may have blades, edges or mechanical devices that can cause severe crush injury to fingers, limbs etc. Take due care when handling and using the product.



Single Operator Only
The product must be operated by a single person only. More than one person operating the product may introduce additional hazards.



Use Face Protection
Wear appropriate full-face protection and take due care as the product or use of the product may present face and eye hazards.



Use Foot Protection
Wear appropriate foot protection and take due care as the product or use of the product may present foot hazards.



Use Eye / Ear / Head Protection
Wear appropriate eye and / or ear and / or head protection and take due care as the product or use of the product may present eye, hearing, and head hazards.



Running Hazard
Do not run on or near the product as doing so may present a fall hazard.



Diving Hazard
Do not dive into the product as doing so may present a neck / head injury hazard.



Adult Supervision Required
Always supervise children and other users of a product to prevent drowning or injury.



Skin Penetration / Puncture Hazard
The product may produce pressure, emit liquids or objects that can cause severe injury to fingers, limbs, blood etc. Take due care when handling and using the product.



Hot Surface Hazard
Be aware that the product may produce high temperatures and hot surfaces that can cause burn injuries.



Flying Debris Hazard
Be aware that the product or use of the product may present hazards produced by flying debris. Wear appropriate clothing and protective devices.



Moving Parts Hazard
Be aware that the product contains or uses mechanical devices that move or rotate. Always wait for moving parts to stop fully before handling the product, adjusting, maintenance etc.

© Equipmed



Carbon-Monoxide Hazard
Do not use the product in confined areas or without adequate ventilation. Carbon-monoxide poisoning can be fatal.



Pull Hazard
Be aware that the product contains or uses mechanical devices that can pull in objects and can cause severe injury to fingers, limbs etc. Take due care when handling and using the product.



Slope / Fall Injury Hazard
Be aware that using the product on sloping surfaces or in slippery conditions may present additional dangers from falls and contact with blades, moving parts, hot surfaces etc.



"Slam Dunk" Warning
Do NOT attempt "slam dunk" manoeuvres as this may result in severe injury due to falling, product breakage or collapse etc.



Electrocution / Electrical Shock Hazard - Outdoor
High voltage or high current electricity may be present or required by the product. Do NOT use in rain, damp, or wet conditions. Electrical shock can be fatal.



Electrocution / Electrical Shock Hazard - Disconnect
High voltage or high current electricity may be present or required by the product. Always disconnect the product from the electrical supply before handling the product, adjusting, maintenance etc.



Power Line Electrocution Hazard
High voltage / high current power lines may be present. Use extreme caution to avoid contact or interference with power lines. Electrical shock can be fatal.



"Kick-Back" Hazard
High level of "kick-back" hazard that can cause the machine to suddenly rotate towards operator. Kick-back injury can be fatal.



Winch Operator Position Hazard
Do NOT stand between winch and load. Do NOT use winch to move people.



Winch Lift Hazard
Do NOT LIFT load vertically. Use machine to PULL only.



Cable Hazard
Ensure that load bearing cable is not kinked or knotted.



Winch Cable Hazard
Ensure that there is a minimum number of cable coils on winching mechanism.



Winch Hook Hazard
Carry hook to load – do NOT throw or run.



Flash / Blinding Hazard
Wear appropriate eye protection for welding. Direct exposure to weld arcs may cause permanent eye injury.



Laser Hazard
Laser may be in use – do NOT look directly at laser or allow others to.



Choking Hazard
Children under 8 yrs. Can choke or suffocate on small parts. Adult supervision required.

© Equipmed


Table of Contents

Safety **2**
Safety Symbols 6
Installation **9**
Operation **11**
Charging Process 11
Hand Control Features 11
Hand Control Operation 11
Starting Use 12
Storage 12
Carrying 13
Emergency Stop 13
Environmental Impact on Equipment 13
Lithium Battery Care and Maintenance **14**
Top 3 Battery Care Requirements 14
Maintenance **16**
After Use and Storage 16
Troubleshooting **17**
Specifications **18**


Installation

Before use, check all parts for damage or defects. Do not use the bath lift if any parts are damaged or missing.


1. After unpacking the bath lift, read this manual before installation.
2. Remove the backrest and seat base from the packaging.
3. Install the backrest onto the seat base as shown.
4. Align the backrest hook with the base guide bar.

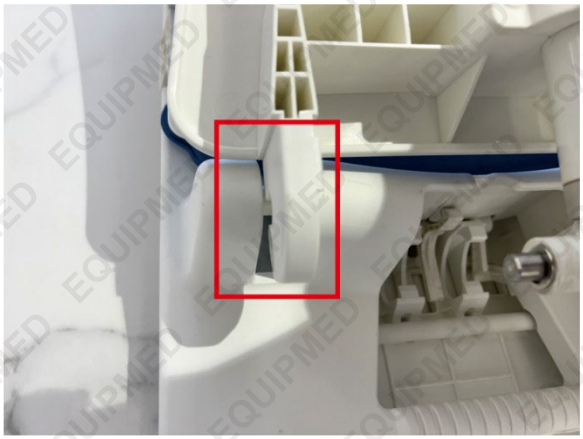

5. Press the hooks on each side of the backrest to secure them to the base.

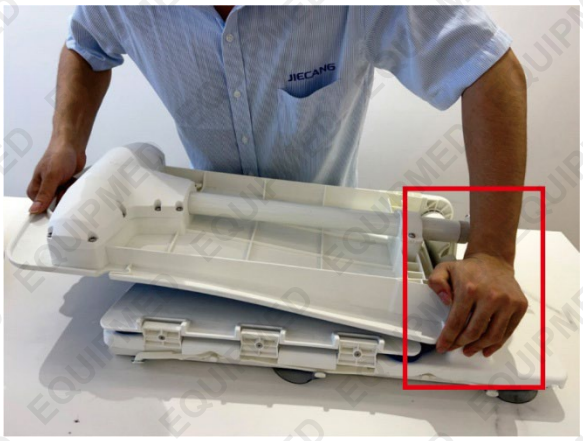

6. Check that the backrest and base are securely connected.

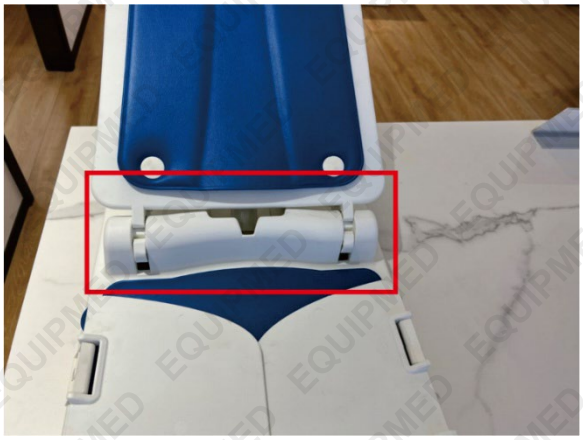

7. Ensure the top clevis of the actuator is correctly connected to the base hook lock.



CORRECT √



WRONG ×

8. Check that the six suction cups on the underside of the seat have not loosened during transport or handling.
9. If any suction cups are loose, tighten them before use.
10. Check that there are no obstructions in the bath that may cause the bath lift to tip over or affect its lifting function, including:
    - shower heads
    - taps
    - fittings
    - small loose items
11. Place the bath lift into the bath. Ask another person for help if required.
12. Ensure the floor is dry and free from obstructions before handling or placing the bath lift into the bath.
13. When placing the bath lift into the bath:
    - ensure there is enough space behind the bath lift
    - prevent the bath lift from hitting the bath or other bathroom objects
    - ensure taps, shower heads and other fittings are outside the lift's operating range
14. Ensure all six suction cups are in full contact with the bottom surface of the bath.
15. Before connecting the hand control, ensure the bath lift plug and hand control plug are completely dry.
16. Position the hand control cable outside the bath lift to prevent it from being pressed, pulled or damaged.
17. Before each use, check that the bath lift raises and lowers properly and safely.

Operation

Charging Process

1. Disconnect the hand control from the bath lift.
2. Ensure the connectors and charger are completely dry.
3. Connect the female plug of the hand control to the male plug of the charger.
4. Plug the charger into a suitable power point.
5. Charging voltage range: 100–240 V AC.
6. When charging is complete, disconnect the charger promptly to help maintain battery life.

Hand Control Features



1. **Upper triangle button**: raises the bath lift.
2. **Lower triangle button**: lowers the bath lift.
3. **Round button**: emergency stop button. Press this button to stop the bath lift immediately. If the bath lift stops while lowering, press the up button first, then press the down button to lower it again.
4. **Indicator light**: shows the operating status of the bath lift.

Hand Control Operation

1. Press the up button to raise the bath lift.
2. Press the down button to lower the bath lift.
3. Press the stop button to activate the emergency stop function.
4. If the stop button is pressed while the bath lift is lowering, the bath lift will stop and will not continue lowering until the down button is pressed again.
5. The hand control has textured buttons to assist users with visual impairments.
6. The hand control operates in press-and-hold mode:
   - the button works while it is pressed
   - the button stops working when it is released

© Equipmed

Starting Use

1. Open the back of the bath lift and connect the hand control plug to the socket at the back of the unit.






2. Raise the bath lift until the side flaps rest securely on both sides of the bath.
3. With your back facing the bath, sit on the outside side flap.
4. While seated, move carefully to the centre of the bath lift seat.
5. Place both feet into the bath.
6. Once seated comfortably in the bath, press and hold the down button to lower the bath lift to the desired position.
7. Keep hands and feet away from moving parts during operation.
8. If the down button is held continuously, the bath lift will stop automatically at its lowest position.
9. Where possible, the user should operate the hand control instead of the assistant or carer. This may help the user feel more comfortable and reduce risk during use.

Storage

- After use, the user should leave the bath and wait for the water to drain before lowering the bath lift to its lowest position and folding the backrest.
- Before lifting the bath lift, release the suction cups by gently lifting the side flaps or armrests.

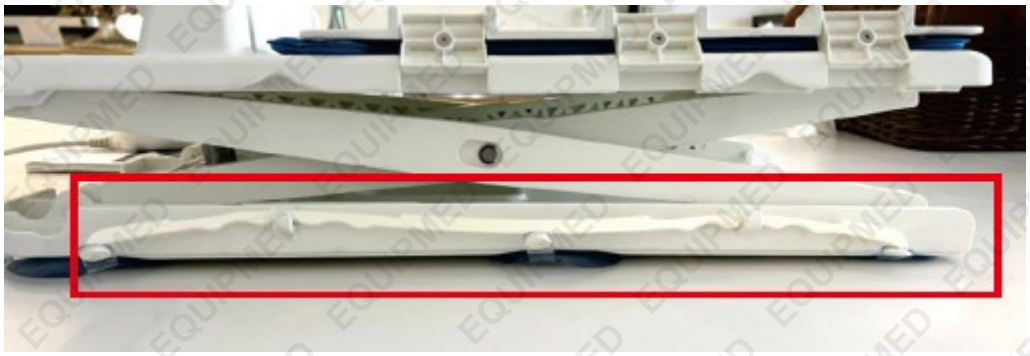

- Rinse the bath lift with clean water to remove foam, stains and residue.

- Keep the hand control connected to the bath lift plug while rinsing the bath lift.
- Disconnect the hand control from the bath lift and wipe it clean with a dry cloth.
- Store the bath lift and hand control in a clean, dry place.

Carrying

- Where possible, carry the backrest and base separately.





- Carry the whole bath lift only when it is safe to do so.
- Use correct lifting techniques to avoid injury.

Emergency Stop

- Press the stop button if the bath lift does not stop normally.
- If the emergency stop button is pressed while the bath lift is lowering:
  - the down button will be locked
  - press the up button first
  - press the down button again to continue lowering the bath lift

Environmental Impact on Equipment

- Keep the bath lift away from direct sunlight to help prevent material ageing.
- Keep the bath lift and battery away from naked flames and high heat to help prevent battery damage, material ageing and deformation.
- Keep the charger away from high-humidity areas to reduce the risk of short circuiting.

© Equipmed

Lithium Battery Care and Maintenance

Top 3 Battery Care Requirements

1. **Charge battery before first use**. The battery may arrive with some charge, regardless of this the battery must be charged for 2–3 hours on the initial charge and 2–3 hours for each subsequent charge.



2. **Charge battery immediately when it is depleted.** Do not store a depleted battery. This will permanently impact the life of the battery. If the battery is going to be stored for a period of time, fully charge the battery before storing.



To ensure the battery maintains optimum performance, follow the guidelines outlined below.

- For best results, charge batteries in temperatures between 10°C and 30°C
- When battery is charged, disconnect the charger from the power supply and remove the battery from the charger.
- If a battery is stored or otherwise unused for an extended period, and the battery has no charge remaining when you check it, consider it to be damaged. Do not attempt to recharge it or to use it. Replace it with a new battery.
- Use only with the battery charger specified by the manufacturer.
- Ensure that the battery charger electrical cord is not subjected to damage or stress. Do not operate the battery charger if it has a damaged electrical cord or plug.
- Do not use the charger in wet areas or expose it to rain or water.
- The typical estimated life of a Lithium-Ion battery is about two to three years or 300 to 500 charge cycles, whichever occurs first. One charge cycle is a period of use from fully charged, to fully discharged, and fully recharged again.
- Explosion hazard – protect the battery against heat, for example, direct sunlight and fire. Do not store batteries in vehicles or locations subject to heat.
- Batteries can explode in the presence of an ignition source, such as a pilot light. To reduce the risk of serious personal injury, never use any cordless product in the presence of open flame. Exploding batteries can propel debris and chemicals.
- Do not open the battery – danger of short-circuiting and/or explosion.
- Do not attempt to destroy or disassemble the battery pack or remove any of its components.
- Do not touch the battery terminals with metal objects and/or body parts as short-circuit and/or personal injury may result.
- Explosion hazard – do not open and/or short-circuit the battery.
- Poison hazard – battery leakage (liquid ejection). Under extreme conditions, liquid may be ejected from the battery – avoid contact. If contact accidentally occurs, immediately seek medical attention, and flush with water. If liquid contacts eyes, immediately seek medical attention, flush with clean water for at least 10 minutes. Liquid ejected from the battery may cause irritation or burns.
- Dispose of non-serviceable batteries in an environmentally responsible manner and according to local regulations.

Failure to adhere to the requirements above may result in the battery becoming inoperative or battery life being reduced. In extreme cases or abuse, battery overheating, or fire may result.

© Equipmed

DANGER - BATTERY CHARGING, USE AND STORAGE

Like mobile phones, laptops, cordless vacuum cleaners etc, this product may contain or be related to a lithium ion battery pack (lithium ion batteries and lithium ion battery related items only). Lithium ion batteries allow fast charging and high capacities, however care must be taken to reduce chances of a fault, and increase safety in case a fault occurs. Catastrophic incidents involving lithium ion batteries have a tendency to result in large amounts of smoke and fire. Be sure to follow the points below to reduce the likelihood of danger whenever you use, charge or store any lithium ion battery.

- Never charge a battery that has been damaged or subjected to a forceful impact.
- Never charge a battery that has been exposed to moisture or water.
- Never charge a battery with a charger that wasn't supplied with the battery.
- Never use a battery in a device that it wasn't specifically designed for.
- Never use a battery with a modified or damaged device.
- Never modify the battery, product, charger, circuitry or wiring.
- Never charge a battery in an exit path that may be used by people to leave a room or building in case of an emergency, for example, the hallway or foyer of a residence.
- Never charge a battery in an area or building that doesn't have the recommended number of working smoke detectors that have been tested in accordance with local guidelines.
- Never charge a battery when the battery or charger is displaying a problem or charging incorrectly.
- Always charge batteries on a non-combustible surface, for example, on a concrete floor.
- Always charge batteries at least 3m away from anything combustible, such as floors, wall surfaces, furniture, storage shelves etc.
- Always charge batteries in an open area (outside if possible) ensuring the battery and charger are away from moisture and water.
- Always disconnect batteries from the charger after charging.
- Always remove or disconnect batteries from the device/product when not in use.
- Always ensure batteries are charged under adult supervision. Do not charge unattended. Do not charge batteries while no one is home or while people are sleeping.
- Keep batteries away from heat and moisture.
- Damaged batteries should be disposed of through your local e-waste collection service.
- Avoid mixing and matching chargers. Unplug products when fully charged, and charge batteries in a cool, dry place away from combustible materials like beds, lounges, or carpet.
- During charging and use, check your lithium-ion batteries for overheating, signs of swelling, leaking, or venting gas. Immediately stop using your product if these signs are present.
- Allow batteries to cool after use and before charging.

Maintenance

- After each bath, rinse the bath lift with warm water.
- If the bath lift is dirty, clean it with a cloth or soft brush, then dry it with a cloth.
- Do not clean the hand control or actuator with detergent or disinfectant.
- If necessary, wipe the hand control and actuator with a damp cloth using a small amount of water.
- Do not immerse the hand control or actuator in water.
- Do not use high-concentration disinfectant to clean the bath lift.
- When cleaning the lower part of the seat, take care around the scissor mechanism to avoid pinching or trapping your hands.
- Clean the bath lift thoroughly between different users.
- Ensure the bath lift is completely dry before storage.
- Ask another person to assist with cleaning and maintenance if required.
- Unplug the charger from the power point before cleaning.
- Do not rinse the charger with water. Clean it only with a dry cloth.
- After each use:
  - wipe off water droplets with a dry cloth
  - allow the bath lift to dry in a well-ventilated area
  - store the bath lift in its box or in a clean, dry place
- Every 3 months:
  - wipe the bath lift with 75% medical alcohol
  - disinfect all external surfaces
  - allow the bath lift to dry in a well-ventilated area for 1 hour
  - store the bath lift in its box or in a clean, dry place
- The bath lift must be inspected and serviced by a qualified professional every 12 months.

After Use and Storage

- If the bath lift has been stored in a cold environment, allow it to warm naturally to room temperature before use.
- After each bath, rinse the bath lift with clean water and allow it to dry naturally.

Troubleshooting

If the bath lift does not raise or lower when unloaded:

- listen for motor noise
- if the motor runs but the bath lift does not move, the actuator gear may be worn or damaged
- if there is no motor noise, the internal drive system may be blocked
- check whether water has entered the unit
- if the issue continues, contact your dealer or the manufacturer

If the bath lift is unstable in the bath:

- check that all suction cups are firmly attached to the bottom of the bath
- ensure the suction cups are clean and in full contact with the bath surface

If the bath lift operates normally when empty but will not raise with the user seated:

- check that the user does not exceed the maximum load limit
- the bath lift may have entered overcurrent protection
- remove the load and check the bath lift before using it again

If the bath lift does not raise or lower when empty:

- check whether the buzzer alarm sounds
- if the buzzer sounds, the battery power is low and the hand control must be charged
- if there is no response, contact your dealer or the manufacturer

Specifications

Max User Load	136kg
Seat Material	PVC
Seat Size	42 x 69cm
Recline Angle	90°–140°
Remote Control	Wired remote control
IP Rating	IP68
Battery	Rechargeable 2200mAh Lithium Coin Cell
Charge Time	2–3 hours
Non-Slip Base	Silicone and plastic
Assembly Required	Yes
Lift Mechanism	Electric Worm Drive
Lift Action	Scissor

© Equipmed

Some experts believe that the incorrect or prolonged use of almost any product may cause serious injury or death. To help reduce your risk of serious injury or death, refer to the information below. For more information, see www.datastreamserver.com/safety

- Consult all documentation, packaging, and product labelling before use. Note that some products feature documentation available online. It is recommended to print and retain the documentation.
- Before each use, check the product for loose/broken/damaged/missing parts, wear, or leaks (if applicable). Never use a product with loose/broken/damaged/missing parts, wear, or leaks.
- Products must be inspected and serviced (if applicable) by a qualified technician every 6 months. This is based on average residential use by persons of average size and strength, and on a property of average metropolitan size. Use beyond these recommendations may require more frequent inspections/servicing.
- Ensure that all users of the product have completed a suitable industry recognised training course before being allowed access to the product.
- The product has been supplied by a general merchandise retailer that may not be familiar with your specific application or description of application. Be sure to attain third-party approval from a qualified specialist for your application before use, regardless of any assurances from the retailer or its representatives.
- This product is not intended for use where fail-safe operation is required. As with any product (for example, automobile, computer, toaster), there is the possibility of technical issues that may require the repair or replacement of parts, or the product itself. If the possibility of such failure and the associated time it may take to rectify could in any way inconvenience the user, business, or employee, or financially affect the user, business, or employee, then the product is not suitable for your requirements. This product is not intended for use where incorrect operation or a failure of any kind, including but not limited to, a condition requiring product return, replacement, parts replacement, or service by a technician may cause financial loss, loss of employee time or an inconvenience requiring compensation.
- If this product has been purchased in error when considering the information presented here, contact the retailer directly for details of their returns policy, if required.

© Equipmed. All rights reserved. No part of this document, including descriptive content, concepts, ideas, diagrams, or images may be reproduced or transmitted in any form or by any means, electronic or mechanical, including photocopying, scanning, or recording, or any information storage and retrieval system, without express permission or consent from the publisher.